

The Great White® Ultra Bur Ideal for Crown Preparation

Smooth tip helps ensure marginal integrity

Dentated carbide specifically designed for superb cutting action

Cuts fast with minimal clogging

Complete the bulk reduction process with one instrument

GREAT WHITE®
ULTRA
CROWN & BRIDGE PREPARATION CARBIDE BURS

"The Great White® Ultra is a versatile bur designed for crown and bridge. It cuts smoothly and efficiently through dentin, enamel, amalgam, composite and cast metal restorations. The unique blade geometry ensures smooth cutting action even when moving from dentin to amalgam or cast metal restorations."


- Steven Abrams, DDS

Traditional diamond instruments often clog during the bulk reduction task in crown preparation. The diamond particles tend to clog with debris, particularly when cutting non-brittle substrates like amalgam. There is no clear way to wash away the debris. As a result diamonds quickly lose their cutting efficiency.

Unlike a diamond, the Great White® Ultra bur provides a faster cutting speed, with less operator pressure resulting in easier instrument advancement.


The Great White® Ultra Bur vs. Diamond Instruments

- Cuts faster (up to 40%) and with significantly less effort than diamonds
- Cuts restorative material faster and easier and transitions between different substrates with minimal clogging and no drag
- Results in less chair time and greater efficiency for the practice
- Bladed cutting tip leaves smooth precise margins
- Less operator fatigue due to feather touch application
- Save up to 4 minutes when cutting a 10 minute crown preparation*

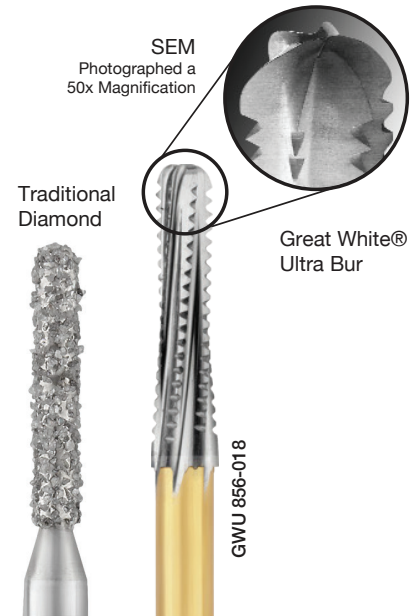


Versatile Performance

Our 856 Series is designed to create a rounded axial-gingival line angle for metal to ceramic restorations.

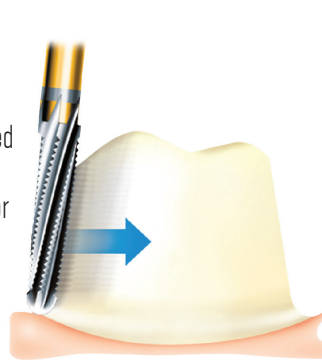

Our 847 Series is designed to develop a 90° axial-gingival line angle suitable for anterior all ceramic and PFM restorations.

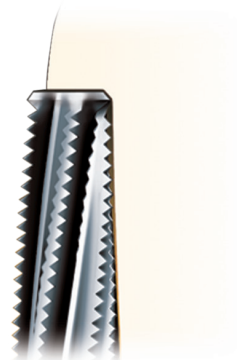

100% Satisfaction Guaranteed

*Market research and test studies available upon request. ©2022 SS White All Rights Reserved. SS White® and Great White® Ultra are registered trademarks of SS White. †SS White's 100% guarantee is on all of our products. If you are not completely satisfied with your purchase, we will issue a full refund of the purchase price through your distributor.

SS WHITE Dental®

Better Patient Outcomes
Improved Efficiency
Faster Practice Growth

1145 Towbin Avenue Lakewood, New Jersey 08701

<https://www.sswhitedental.com/gw-samples>

GREAT WHITE[®] ULTRA

CROWN & BRIDGE PREPARATION CARBIDE BURS

FOOTBALL



GWU 379-023
18193-5

5 PK
HEAD LENGTH (mm) 4.2
MAX HEAD DIAMETER (mm) 2.3

SHORT FLAT END TAPER



GWU 845KR-016 18195-5
GWU 845KR-018 18196-5
GWU 845KR-025 18197-5

5 PK
HEAD LENGTH (mm) 4.0
MAX HEAD DIAMETER (mm) 1.6

FLAT END TAPER



GWU 847-016 18200-5
GWU 847-018 18191-5
GWU 847-020 18192-5

5 PK
HEAD LENGTH (mm) 8.0
MAX HEAD DIAMETER (mm) 1.6

FLAT END TAPER SHORT SHANK



GWU 847-016S 18186-5
GWU 847-018S 18187-5
GWU 847-020S 18185-5

5 PK
HEAD LENGTH (mm) 8.0
MAX HEAD DIAMETER (mm) 1.6

ROUND END TAPER



GWU 855-025 18194-5
GWU 856-016 18198-5
GWU 856-018 18188-5
GWU 856-020 18189-5
GWU856-023 18190-5

5 PK
HEAD LENGTH (mm) 7.0
MAX HEAD DIAMETER (mm) 2.5

ROUND END TAPER SHORT SHANK




GWU 856-016S 18182-5
GWU 856-018S 18183-5
GWU 856-020S 18184-5

5 PK
HEAD LENGTH (mm) 8.0
MAX HEAD DIAMETER (mm) 1.6

Crown and Bridge Preparation Kit

CAD/CAM Preparation & Finishing Kit

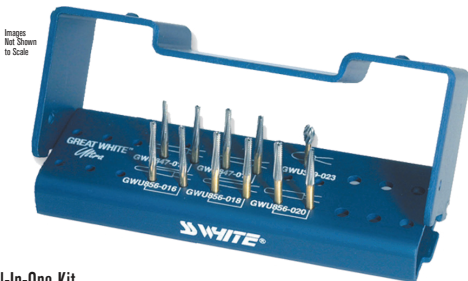

All-In-One Kit
Designed for the efficient atraumatic removal of failed amalgams, composites, zirconia and PFM Crowns. **Kit Reorder No. 18150**

Kit Contents:
2-Great White[®] Ultra #856-016 - **#18198-5**
2-Great White[®] Ultra #856-018 - **#18188-5**
2-Great White[®] Ultra #856-020 - **#18189-5**
1-Great White[®] Ultra #379-023 - **#18193-5**
2-Great White[®] Ultra #847-016 - **#18200-5**
2-Great White[®] Ultra #847-018 - **#18191-5**
1-Autoclavable Bur Block

The Great White[®] Ultra:

- Cuts quickly and smoothly through enamel
- Negotiates amalgam and other restorative materials with minimal clogging and no "drag or stalling"
- Allows completion of the bulk reduction task in the crown preparation with a single instrument
- More efficient –less chair time
- Helps provide precise margins
(Please reference SEM photograph)



Add efficiency and organization to CAD/CAM preparations with a streamlined selection of rotary instruments for milled ceramic crowns, inlays and partial coverage restorations. Featuring the innovative SS White[®] Great White[®] 1.5mm & 2mm depth cutter along with selected carbide burs, diamond instruments and Jazz[®] Polishers, organized in an autoclavable metal bur block. **Kit Reorder No. 64502**

Kit Contents:
Great White Depth Cutter 1 **#15040**
Great White Depth Cutter 1 **#15041**
FG #169 **#16884**
Finishing Bur #7675 - **#15675-5**
Great White[®] Ultra #379-023 - **#18193-5**
Great White[®] Ultra #856-016S - **#18182-5**
Piranha Diamond #881-016 - **#881-016VF**
Jazz[®] Polisher Fine Flame P3S - **#89012**
Jazz[®] Polisher Med Flame P3S - **#89008**
Jazz[®] Polisher Coarse Flame P3S - **#89004**

Great White[®] Depth Cutters



GWDC1 15040
STERILE ROLL 25 PK 10933
HEAD LENGTH (mm) 1.5mm

GWDC2 15041
10138
2mm

Great White Depth Cutters are specifically designed with a self-limiting stop to produce 1.5 mm or 2 mm occlusal depth cuts necessary for CAD/CAM crown, inlay or onlay preparations. Great White[®] Depth Cutters will not clog when used on amalgam and composite restorations, compared to diamond depth cutters, making it an ideal universal instrument when initiating preparations on restored or non-restored occlusal surfaces.