




CRA  
CONFIRMED  
USEFUL

## NOTEWORTHY —

### FISSURE EXCAVATING BURS — SMALLER SIZES PRODUCE MORE CONSERVATIVE CUT

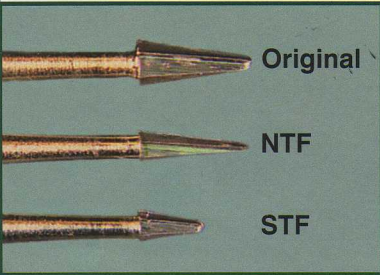

\$3.50 each

SS White Burs, Inc.  
1145 Towbin Ave • Lakewood, NJ 08701 • U.S.A.  
732-905-1100 • 800-535-2877 • Fax: 732-905-0987  
www.sswiteburs.com

#### FISSUROTOMY BURS—

#### NTF (Narrow Tapered Fissure) & STF (Short Tapered Fissure)

6-bladed carbide burs in 3 diameters (1.1mm, 0.7mm, 0.6mm), 2 lengths (2.5mm, 1.5mm), & standard or short shanks. Designed to excavate pits & fissures, & interproximal caries, & remove small old amalgam restorations. Advantages: (1) Bur design can produce very conservative occlusal preps when used properly; (2) Can be used in place of air abrasion; (3) Anesthetic often not required, especially when cutting only enamel; & (4) Sharpness produces fast cut. Disadvantages: (1) Cost & (2) Like all carbide burs, steam sterilization dulls & corrodes. For best results a new bur should be used for each patient.

**83%** of 18 CRA Evaluators stated this product would replace products they use currently, & 100% rated it excellent or good & worthy of trial by colleagues.

